tore at

9879

SimpleChIP® Human DNAJB9 Exon 1 Primers

Cell Signaling
TECHNOLOGY®

Support: +1-978-867-2388 (U.S.) www.cellsignal.com/support

Orders: 877-616-2355 (U.S.) orders@cellsignal.com


Entrez-Gene ID #4189 UniProt ID #Q9UBS3

Small 500 µl (250 PCR reactions)

For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Species Cross-Reactivity	Primer Anneal/Extension	PCR Product Length
ChIP	H	65°C	100 bp


Description: SimpleChIP® Human DNAJB9 Exon 1 Primers contain a mix of forward and reverse PCR primers that are specific to exon 1 of the human DnaJ homolog subfamily B member 9 gene. These primers can be used to amplify DNA that has been isolated using chromatin immunoprecipitations (ChIP). Primers have been optimized for use with SimpleChIP® Universal qPCR Master Mix #88989 and have been tested in conjunction with SimpleChIP® Plus Enzymatic Chromatin IP Kits #9004 and #9005 and ChIP-validated antibodies from Cell Signaling Technology®.

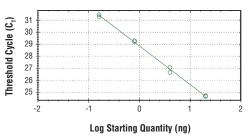

New 10/18

E = 104.6% R² = 0.997

SimpleChIP® Human DNAJB9 Exon 1 Primers were tested on DNA isolated from cross-linked cells using the SimpleChIP® Plus Enzymatic Chromatin IP Kit (Magnetic Beads) #9005. Real-time PCR was performed in duplicate on a serial dilution of 2% total input DNA (20 ng, 4 ng, 0.8 ng, and 0.16 ng) using a real-time PCR detection system and SimpleChIP® Universal qPCR Master Mix #88989. The PCR amplification efficiency (E) and correlation coefficient (R²) were calculated based on the corresponding threshold cycle (C₁) of each dilution sample during 40 cycles of real-time PCR (95°C denaturation for 15 sec, 65°C anneal/extension for 60 sec).

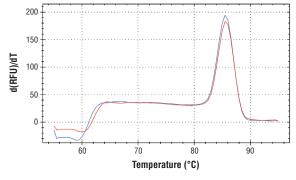

PCR product melting curves were obtained for real-time PCR reactions performed using SimpleChIP® Human DNAJB9 Exon 1 Primers. Data is shown for both duplicate PCR reactions using 20 ng of total DNA. The melt curve consists of 80 melt cycles, starting at 55°C with increments of 0.5°C per cycle. Each peak is formed from the degradation of a single PCR product.

Storage: Supplied in nuclease-free water at a concentration of 5 μ M (each primer is at a final concentration of 5 μ M). *Store at*

Directions for Use:

- 1. Label the appropriate number of PCR tubes or PCR plates compatible with the model of real-time PCR machine to be used. PCR reactions should be performed in duplicate and should include a tube with no DNA to control for contamination, and a serial dilution of a 2% total input chromatin DNA (undiluted, 1:5, 1:25, 1:125), which is used to create a standard curve and determine amplification efficiency.
- 2. Add 2 μ l of the appropriate ChIP DNA sample to each tube or well of the PCR plate.
- Prepare a master PCR reaction mix as described below. Add enough reagents for two extra reactions to account for loss of volume. Add 18 µl of the master PCR reaction mix to each PCR reaction tube or well of the PCR plate.

 $\begin{array}{lll} \textbf{Reagent} & \textbf{Volume for 1 PCR Reaction (20 } \mu \textbf{J}) \\ \textbf{Nuclease-free H}_2 \textbf{O} & \textbf{6} \; \mu \textbf{I} \\ \textbf{5} \; \mu \textbf{M SimpleChIP}^{\otimes} \; \textbf{Primers} & \textbf{2} \; \mu \textbf{I} \\ \textbf{SimpleChIP}^{\otimes} \; \textbf{Universal qPCR Master Mix \#88989} & \textbf{10} \; \mu \textbf{I} \\ \textbf{Output} & \textbf{All positions} & \textbf{All positions} & \textbf{All positions} \\ \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{All positions} \\ \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{All positions} \\ \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{Output} \\ \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{Output} \\ \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{Output} \\ \textbf{Output} & \textbf{Output} & \textbf{Output} \\ \textbf{Output} & \textbf{Output} & \textbf{Output} & \textbf{Output} \\ \textbf{Output} & \textbf{Output} \\ \textbf{Output} & \textbf{Output} \\ \textbf{Output} & \textbf{Output} \\ \textbf{$

- 4. Start the following PCR reaction program:
 - a. Initial Denaturation: 95°C for 3 min.
 - b. Denaturation: 95°C for 15 sec.
 - c. Anneal and Extension: Primer-specific temp. for 60 sec.
 - d. Repeat steps b and c for a total of 40 cycles.
- 5. Analyze quantitative PCR results using software provided with the real-time PCR machine.

Thank you for your recent purchase. If you would like to provide a review visit cellsignal.com/comments.

www.cellsignal.com